

TEAC UR-4MD

High-Definition video recorder



EN 60601-1

The UR-4MD is the latest high-definition medical video recorder from TEAC. Designed for a quick set up and ease of use the high-definition recorder offer a lot of features. For example the auto sensing video connection and automatic video resolution detection. The front control panel allows users to intuitively record videos and photos in cristal clear HD. The internal hard disk reliable save recording and data back-ups. In addition, the user can also simultaneously record to USB memory and USB hard drives to archive, share and play back on other devices and locations. Another plus is the networking capability of the UR-4MD.


Benefits

- › Full-HD 1920 x 1080 video recording
- › Extreme ease of use
- › Auto-sensing I/O and video resolution
- › Simultaneous recording to multiple media (up to 2 external USB- drives and 1 internal HD)
- › 3.5" LDC-monitor for Live View, Playback or GUI
- › Networking (FTP and CIFS) RJ45 10/100/1000 Ethernet
- › Auto-shutdown in case of power interruption
- › S.M.A.R.T. monitors/ reports - Self-Monitoring, Analysis and Reporting Technology
- › Fluid resistant front Panel



General	UR-4MD	
Item No.	2011656	
Recordable Devices	USB-Flash-Drive USB- Hard Drive Internal Hard Drive (500 GB)	
Display	3.5 Zoll LCD	
Power supply	100 V-240 V, 50/60 Hz	
Temperature operating	5 to 40 °C	
Temperature storage	-20 to 60 °C	
Humidity operating	5 to 80% (no condensation)	
Dimensions (WHD)	210x 88.5x 235 mm	
Weight	approx 2.6 kg	
Standards	IEC60601-1, CSA60601-1, EN60601-1, CE, IEC60601-1-2, EN60601-1-2, ICES-003, VCCI Klasse B, FCC Klasse B, FDA Klasse 1, ISO13485, CSA NRTL	
Accessories (optional)	Remote control unit, Foot-switch, keyboard, Barcode / Card Reader	
Recording Format		
Still image	Jpeg	
Video	MPEG-4 AVC/H.264	
Input/ Output		
Video - IN/OUT	SDI BNC x 1	High-Definition, NTSC, PAL
	DVI (Single Link) x 1	High-Definition
	S-Video x 1	NTSC, PAL
	Composite (BNC) x 1	NTSC, PAL
Audio - IN/OUT	SD/HD-SDI Thru x 1 (only Video-OUT)	High-Definition, NTSC, PAL
	SDI (BNC) x 1	* SDI-Video input terminal use in common
	Stereo-mini jack x 2	LINE, MIC (only Audio-IN)
Other Interfaces	RJ-45 (10BASE-T / 100BASE-TX / 1000Base-T) x 1	
	USB 2.0 Typ A x 3	2 front connector / 1 rear connector
	USB 2.0 Typ B x 1	1 rear connector for accessing the internal Hard Drive (on UR-4MD)
	Stereo mini jack x 3	for external trigger

System configuration

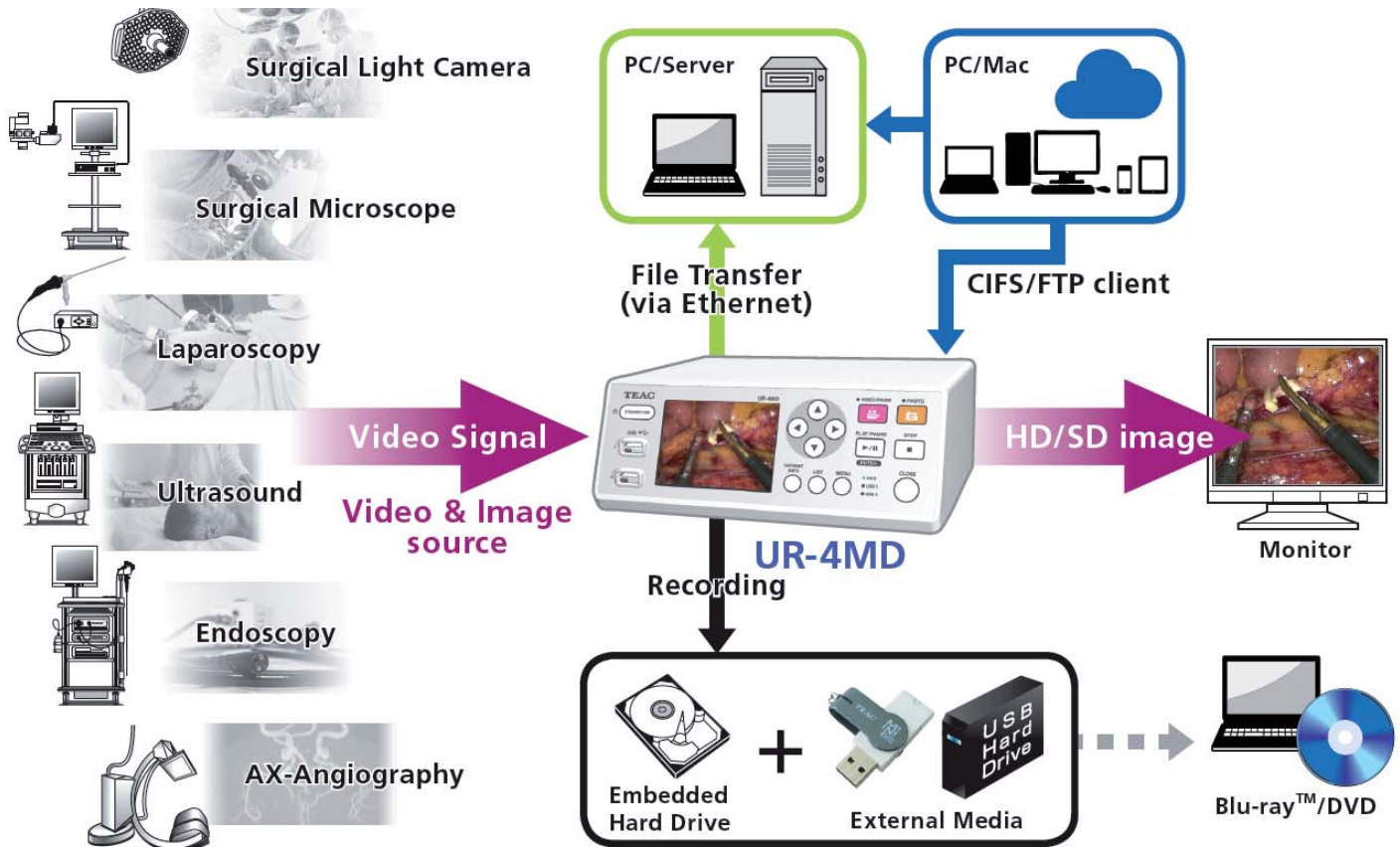

Approximate Recording Times for 500 GB Hard Drive

High-Definition 1920 x 1080

12 Mb/s	18 Mb/s	24 Mb/s
90 hours	60 hours	45 hours

Standard-Definition 1920 x 1080

5 Mb/s	8 Mb/s	10 Mb/s
180 hours	135 hours	100 hours